



LC Docking station for pumps, valves and detectors

With more than 4000 module combinations you will find the perfect configuration for your LC. The docking station can be customized with up to three compact devices. Choose the properties you want to get your perfect Assistant ASM 2.2L for chromatography and beyond.

Pumps

Choose from 15 different pumps with 10 or 50 ml pump heads and with or without pressure sensor. The material used is stainless steel, ceramic or Hastelloy C (for pumps without pressure sensors only).

Valve drive

The universal valve drive identifies valves via RFID technology and enables to read GLP data. All V 4.1 valves, independent of number of ports and position, are supported.

UV detectors

The compact single wavelength UV detector is available in a basic and fibre optics version. The wavelength can be set between 190 - 500 nm.

Freely combine pumps, valves and detectors in one housing




[Contact us](#)




Eluent delivery



Detection

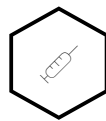

Sample injection

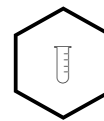

Fraction collection



Sample, eluent &
column selection

Combine three plug-in LC modules in one housing

Liquid Chromatography D...



The plug-in modules can be easily removed by loosening four screws allowing the user to exchange modules in case of service within minutes. Likewise, the configuration of the LC system can be adapted to new requirements. Routine maintenance work e.g. replacing the lamp of a detector are easily performed by the user.

Depending on the integrated modules the ASM 2.2L fulfills many different tasks like eluent delivery, detection, sample and solvent selection, sample injection, column switching or fraction collection. An assistant including a pump, injection valve, and detector features a complete, compact

chromatographic system. As a part of a larger system, the ASM 2.2L is extremely versatile in analytical, preparative and continuous liquid chromatography.

Configurator

AVAILABLE PLUG-IN MODULES



Valve drive

Universal valve drive with
RFID-technology
Controls all V 4.1 valves



UV detectors

Variable single
wavelength UV
detector
190 - 500 nm



Pumps

Compact pumps with or
without pressure sensor
10 or 50 ml pump head

10 or 30 mm pump head
Stainless steel, ceramic,
Hastelloy C

› **See valve drive**

› **See detectors**

› **See pumps**

Configure your LC Assistant

How to find your desired plug-in module combination?

AZURA Assistant ASM 2.2L



The Assistant ASM 2.2L can be equipped with three plug-in modules. Use the assistant configurator to find your desired module combination. Select your plug-in modules for the left, middle and right position in the assistant and get your perfect assistant for chromatography and beyond.

Configuration note

Assistants with following configuration are not allowed:

1) more than two pump modules - a high-pressure gradient is

not support

2) more than one UV detector

3) without plug-in modules

Configurator

That's what the users say



" SMB chromatography is a HPLC technique for the separation of binary mixtures with high productivity and purity. The LC device dock is an appropriate footprint to favorably arrange the pumps and valves in accordance with the SMB process. Second advantage: SMB chromatography stands for a continuous process. The concept of plug-in modules guarantees short maintenance breaks. Since only a small part of the complex capillary construction has to be removed for the replacement of a module in case of service, the downtime of the system is reduced. An obvious benefit is the low space required for the high number of integrated devices."

Dr. Paul Pietsch, SMB Specialist

VALVE DRIVE



Universal valve drive for ASM 2.2L

Assistant module VU 4.1 for valves V 4.1

> Read more

UV DETECTORS



Compact variable single wavelength UV detector for ASM 2.2L



Compact variable single wavelength UV detector with fiber optics for ASM 2.2L

Assistant module UVD 2.1S without flow cell

Assistant module UVD 2.1S fiber optics version without flow cell

> **Read more**

> **Read more**

PUMPS WITH PRESSURE SENSOR



Compact



Compact



Compact



Compact



Compact



Co

<p>pump with pressure sensor, 10 ml/min stainless steel pump head for ASM 2.2L</p> <p>Assistant module P 4.1S, 10 ml, stainless steel</p> <p>> Read more</p>	<p>pump with pressure sensor, 10 ml/min ceramic pump head for ASM 2.2L</p> <p>Assistant module P 4.1S, 10 ml, ceramic</p> <p>> Read more</p>	<p>pump with pressure sensor, 10 ml/min stainless steel pump head for normal phase for ASM 2.2L</p> <p>Assistant module P 4.1S, 10 ml, stainless steel, normal phase</p> <p>> Read more</p>	<p>pump with pressure sensor, 50 ml/min stainless steel pump head for ASM 2.2L</p> <p>Assistant module P 4.1S, 50 ml, stainless steel</p> <p>> Read more</p>	<p>pump with pressure sensor, 50 ml/min ceramic pump head for ASM 2.2L</p> <p>Assistant module P 4.1S, 50 ml, ceramic</p> <p>> Read more</p>	<p>pump with pressure sensor, 50 ml/min stainless steel pump head for normal phase for ASM 2.2L</p> <p>Assistant module P 4.1S, 50 ml, stainless steel, normal phase</p> <p>> Read more</p>
-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------	---------------------------------------------------------------------------------------------------------------------------------------------------------------	--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------	-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------	---------------------------------------------------------------------------------------------------------------------------------------------------------------	--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------


PUMPS WITHOUT PRESSURE SENSOR



**Compact
pump
with 10
ml/min
stainless
steel
pump
head for
ASM
2.2L**



**Compact
pump
with 10
ml/min
ceramic
pump
head for
ASM
2.2L**

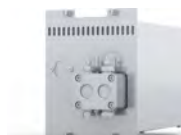


**Compact
pump with
10 ml/min
Hastelloy®
C pump
head for
ASM 2.2L**

Assistant



**Compact
pump
with 50
ml/min
stainless
steel
pump
head for
ASM
2.2L**

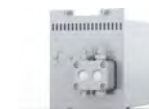

**Compact
pump
with 50
ml/min
ceramic
pump
head for
ASM
2.2L**

Assistant
module P
2.1S, 10
ml,
stainless
steel

> **Read
more**

Assistant
module P
2.1S, 10
ml,
ceramic

> **Read
more**

module P
2.1S, 10 ml,
Hastelloy®
C

> **Read
more**

Assistant
module P
2.1S, 50
ml,
stainless
steel

> **Read
more**

Assistant
module P
2.1S, 50
ml,
ceramic

> **Read
more**

SLOTS WITH NO PLUG-IN MODULES



Empty module for ASM 2.2L

Cover plate for slots without an assistant module

[➤ Read more](#)



Column holder for ASM 2.2L

Column holder for slots without an plug-in module

[➤ Read more](#)

ACCESSORIES



**Display software for AZURA
devices without data acquisition
incl. tablet**

Mobile Control without data
acquisition including tablet and
mount

Article No.: A9607



**Display software for AZURA
devices with data acquisition incl.
tablet**

Mobile Control Chrom with data
acquisition including tablet and
mount

Article No.: A9608

DETAILED INFORMATION

AZURA Assistant ASM 2.2L Instructions (V6803)

Manual English

AZURA Assistant ASM 2.2L Betriebsanleitung (V6803)

Manual German

AZURA Assistant ASM 2.2L Consignes (V6803)

Manual French

List of supported devices in chromatography data systems EN (V1663)

Device support in KNAUER CDS / KNAUER device support in 3rd party software

Information English

Configure your LC Assistant

How to find your desired plug-in module combination?

The Assistant ASM 2.2L can be equipped with three plug-in modules. Use the assistant configurator to find your desired module combination. Select your plug-in modules for the left, middle and right position in the assistant and get your perfect assistant for chromatography and beyond.

Configuration note

Assistants with following configuration are not allowed:

- 1) more than two pump modules - a high-pressure gradient is not support
- 2) more than one UV detector
- 3) without plug-in modules

Previous version: ASM 2.1 L

Here, you will find all information about the LC Assistant ASM 2.1L.

[Go to ASM 2.1L](#)